



Portable HIFU Beauty Machine Simultaneous Control / Independent Control

Our Product Introduction

Basic Information

- Place of Origin: Guangzhou, China
- Brand Name: other
- Certification: CE
- Model Number: hifu-3
- Minimum Order Quantity: 1 Sets
- Price: Discussible
- Packaging Details: Regular Packaging
- Delivery Time: 5-8 working Days
- Payment Terms: L/C, T/T,
- Supply Ability: 10000 pieces/month

HIFU probe and function introduction

Standard 3.0mm and 4.5mm probes

The instrument can be matched with a total of 8 probes with a depth of 1.5~16mm, and different probes can be selected according to the actual situation
2 as standard (up to 8 optional)



Product Specification

- Name: Portable HIFU Beauty Machine
- Body Type: Vertical
- Scope: 5.0mm-25mm
- Size: Portable
- Input Voltage: AC100V-240V
- Working Mode: Simultaneous / Independent Control
- Energy Type: HIFU (High Intensity Focused Ultrasound)
- Lamp Longevity: 300,000 Shots (IPL)
- Spacing: 1.0-10mm(0.5mm/step)
- Highlight: **Portable HIFU Beauty Machine,
HIFU Beauty Machine Simultaneous Control ,
HIFU Beauty Machine Independent Control**

Product Description

Advanced HIFU Beauty Machine With Focus Depth From 5.0mm-25mm

Product Description:

With its vertical body type, this HIFU Beauty Machine is designed for easy use and comfort. It is suitable for all skin types, making it a versatile solution for men and women looking to combat the signs of aging.

One of the standout features of this HIFU Beauty Machine is its working mode – it can operate in simultaneous or independent control modes, depending on your preferences. This allows for a customized treatment plan that suits your individual needs.

The HIFU Anti-wrinkle Device comes in a pack size of 60*44*46cm, making it compact and easy to store. Its specification includes packing at 64 X 81 X 131 Cm, making it easy to transport and use wherever you go.

Investing in the HIFU Anti-wrinkle Device is a great way to achieve younger-looking skin without the need for invasive procedures. Its advanced HIFU technology and customizable treatment modes make it a reliable and effective solution for anyone looking to enhance their skin's appearance.

Features:

Product Name: HIFU Beauty Machine

Supplier Info: Manufacturer

Spacing: 1.0-10mm (0.5mm/step)

Size: Portable

Scope: 5.0mm-25mm

Pack Size: 60*44*46cm

This HIFU Beauty Machine is a portable device that helps reduce wrinkles and tighten skin, making it a great choice for those looking for a non-surgical, non-invasive solution to improve their appearance. With a spacing range of 1.0-10mm and a scope of 5.0mm-25mm, this HIFU Facial Machine is able to target different areas of the face and body with precision. The pack size of 60*44*46cm makes it easy to store and transport.

Technical Parameters:

Working Mode	Simultaneous / Independent Control
Spacing	1.0-10mm(0.5mm/step)
Specification	Packing , 64 X 81 X 131 Cm
Lamp Longevity	300,000 Shots (IPL)
Length	5.0 - 25mm (1.0mm / Step)
Supplier Info	Manufacturer
Pack Size	60*44*46cm
Energy Output	0.1J-2.5J (0.1J/step)
Energy Type	HIFU (High Intensity Focused Ultrasound)
Material	ABS

This Hifu Facial Machine is a Focused Ultrasound Beauty Device that uses HIFU (High Intensity Focused Ultrasound) energy to tighten and rejuvenate the skin. It has a spacing range of 1.0-10mm with 0.5mm/step precision, and a length range of 5.0-25mm with 1.0mm/step precision. The energy output can be adjusted from 0.1J-2.5J with 0.1J/step precision, and the lamp longevity is up to 300,000 shots (IPL). The pack size is 60*44*46cm, and the machine is made of ABS material. It is manufactured by a trusted supplier.

HOW TO LIFT AND TIGHTEN SKIN WITHOUT SURGERY? 5D HIFU, THAT'S HOW.



5D HIFU non-invasive procedure that lifts the neck, chin and brow, and improves lines and wrinkles on the upper chest. Through the ultrasound imaging with Ultrasound energy is delivered precisely to the place where it will benefit you most—for natural-looking results that improve over time.

NATURAL-LOOKING AND NOTICEABLE

The proven results of 5D HIFU are clear, appearing over two to three months as new collagen works to lift skin on the neck, chin, and brow—as well as smooth skin on the upper chest. Results may improve over the following three to six months, and can last a year or more.





HOW MANY TREATMENTS DOES IT TAKE?

The majority of patients only need one treatment; however, some may benefit from more than one treatment, depending on how much skin laxity they have and their body's own biological response to the ultrasound and the collagen-building process. Follow-up HIFU treatments each year may help maintain results.

WHY WE NEED TO DO VAGINAL TIGHTENING?

REASON

As you age there is a natural loss of collagen and muscle tone. This may be experienced as a decrease in firmness of the vaginal tissue/canal. Childbirth can result in vaginal laxity or scarring, resulting in a larger vaginal opening. Strenuous exercise can lead to a weakening of the pelvic floor support muscles. Naturally decreasing levels of estrogen can lead to thinning and dryness of the vaginal skin. Ultimately this can result in a lack of bladder support and leakage of urine with exercise, coughing, laughing, or sneezing.

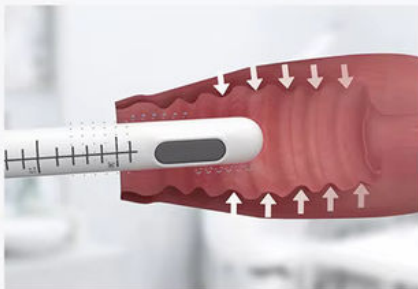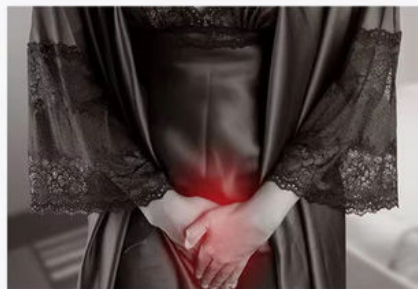

BENEFITS OF 5D HIFU 360° VAGINAL TREATMENT HANDLE:

- Increases Vaginal blood flow and lubrication
- Improves Symptoms of Stress and Urge Incontinence
- Improves pelvic floor strength and tone
- Enhances sexual arousal and orgasm
- Improves Aesthetic Appearance of External Genital Tissue

HIFU probe and function introduction

Standard 3.0mm and 4.5mm probes

The instrument can be matched with a total of 8 probes with a depth of 1.5~16mm, and different probes can be selected according to the actual situation

2 as standard (up to 8 optional)



Product Accessories



Host*1



Probe*1



Standard 3.0mm
4.5mm probe



Power cord*1



4 acrylic sheets
2.5 mm/2 pieces
3.5 mm/2 pieces



6 pieces


Instrument probe function introduction

(Standard 3.0/4.5) (Optional 1.5/6.0/8.0/10.0/13.0/16.0)

Rejuvenate the skin and break the "wrinkle" of the skin

Through high-frequency radio waves, the heat penetrates into the inner layer of the skin, awakening the vitality of collagen and caring for the skin.

III

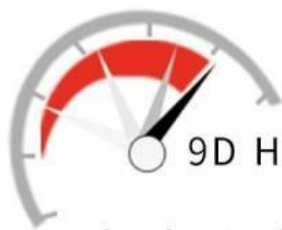

9D HIFU



1 minute / 12 rows of RBI
Multi-line dots save half the time

Up to 12 rows of points at a time



9D HIFU



Hit 12 points at a time

Ordinary HIFU



Hit 1 row at a time

Product parameters



Product name: 9D HIFU	Working mode: pulse digital output mode
Voltage: 110V~220V 60Hz	Rounds: 20000 rounds
Number of probes: 2 as standard (up to 8 optional)	Instrument weight: 14Kg
Energy: 0.2J-2.0J (adjustable 0.1J/sec)	Instrument size: 42*47*30 (cm)
Probe size: 1.5, 3.0, 4.5, 6.0, 8.0, 10.0, 13.0, 16.0 (mm)	

Applications:

The HIFU-3 Facial Lifting Beauty Device is a non-invasive skin tightening beauty device that is designed to provide a range of benefits to users. This HIFU machine is manufactured by other in Guangzhou, China and comes with a CE certification.

The HIFU-3 is a portable device that delivers energy output ranging from 0.1J to 2.5J (0.1J/step). It operates in simultaneous / independent control working mode, and its scope ranges from 5.0mm to 25mm.

With its advanced features and design, the HIFU-3 is ideal for a wide range of product application occasions and scenarios, including:

Skin Lifting and Tightening: The HIFU-3 is perfect for those looking to tighten their skin and reduce the appearance of fine lines and wrinkles. It can be used to lift sagging skin on the forehead, cheeks, chin, and neck.

Body Contouring: This HIFU machine can also be used to contour the body, reducing the appearance of cellulite and tightening skin around the arms, abdomen, and thighs.

Postpartum Recovery: For women who have recently given birth, the HIFU-3 can help restore the body's elasticity and reduce the appearance of stretch marks.

Skin Rejuvenation: The HIFU-3 is also effective in rejuvenating the skin, giving it a more youthful and radiant appearance. It can be used to reduce the appearance of age spots, sun damage, and other skin blemishes.

Beauty Salon and Spa Treatments: The HIFU-3 can be used in beauty salons and spas to provide high-quality skin tightening and body contouring treatments to clients.

The HIFU-3 Facial Lifting Beauty Device is available for purchase with a minimum order quantity of 1 set. The price is discussible, and regular packaging is included. Delivery time is 5-8 working days, and payment terms include L/C and T/T. The supplier is the manufacturer, and the supply ability is 10000 pieces/month.

Overall, the HIFU-3 provides a range of benefits to users and is a versatile and effective tool for those looking to improve the appearance of their skin and contour their bodies.

Customization:

Support and Services:

The HIFU Beauty Machine is a non-invasive cosmetic treatment device that uses high-intensity focused ultrasound technology to lift and tighten the skin. It is designed to improve the appearance of fine lines, wrinkles, and sagging skin on the face and neck.

Our product technical support team is available to assist with any questions or issues you may have with the HIFU Beauty Machine. We offer comprehensive training and support materials to ensure that you can use the device safely and effectively.

In addition to technical support, we also offer a range of services to help you get the most out of your HIFU Beauty Machine. These

include marketing materials and guidance to help you promote and sell the device to your customers, as well as ongoing product updates and enhancements to ensure that your device is always up-to-date and operating at peak performance.

Packing and Shipping:

Product Packaging:

1 HIFU Beauty Machine
1 Power Adapter
1 User Manual
1 Carrying Case

Shipping:

Shipping Method: FedEx
Estimated Delivery Time: 3-5 Business Days
Shipping Cost: \$20
Tracking Number: Provided via Email

FAQ:

Q: What is the brand name of the HIFU Beauty Machine?

A: The brand name of the HIFU Beauty Machine is other.

Q: What is the model number of the HIFU Beauty Machine?

A: The model number of the HIFU Beauty Machine is hifu-3.

Q: Where is the HIFU Beauty Machine manufactured?

A: The HIFU Beauty Machine is manufactured in Guangzhou, China.

Q: Is the HIFU Beauty Machine certified?

A: Yes, the HIFU Beauty Machine is CE certified.

Q: What is the minimum order quantity for the HIFU Beauty Machine?

A: The minimum order quantity for the HIFU Beauty Machine is 1 set.

Q: What are the payment terms for the HIFU Beauty Machine?

A: The payment terms for the HIFU Beauty Machine are L/C and T/T.

Q: What is the supply ability of the HIFU Beauty Machine?

A: The supply ability of the HIFU Beauty Machine is 10,000 pieces per month.

Q: What is the packaging details of the HIFU Beauty Machine?

A: The packaging details of the HIFU Beauty Machine is regular packaging.

Q: What is the delivery time for the HIFU Beauty Machine?

A: The delivery time for the HIFU Beauty Machine is 5-8 working days.

Q: Is the price of the HIFU Beauty Machine negotiable?

A: Yes, the price of the HIFU Beauty Machine is discussible.



Guangzhou Lizhirun Technology Development Co., Ltd



18688905627



cenglinlin25@gmail.com



hifubeaity-machine.com

Room 502, 5th Floor, Building F, Huixin Industrial Zone, No. 89, Hetai Road, Hebain Village, Baiyun District,
Guangzhou