




Ultrasound HIFU Beauty Machine ABS Ultrasound Beauty Device With Adjustable Spacing

Our Product Introduction

For more products please visit us on hifubeauty-machine.com

Basic Information

- Place of Origin: Guangzhou, China
- Brand Name: other
- Certification: CE
- Model Number: hifu-3
- Minimum Order Quantity: 1 Sets
- Price: Discussible
- Packaging Details: Regular Packaging
- Delivery Time: 5-8 working Days
- Payment Terms: L/C, T/T, D/A, D/P, Western Union, MoneyGram
- Supply Ability: 10000 pieces/month

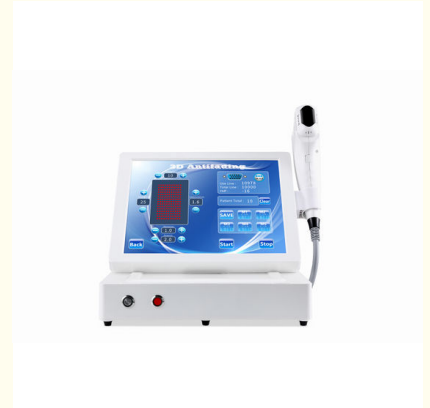

Product Specification

- Energy Type: HIFU (High Intensity Focused Ultrasound)
- Material: ABS
- Size: Portable
- Suitable Skin Types: All
- Spacing: 1.0-10mm(0.5mm/step)
- Input Voltage: AC100V-240V
- Specification: Packing , 64 X 81 X 131 Cm
- Length: 5.0 - 25mm (1.0mm / Step)
- Highlight: **Ultrasound HIFU Beauty Machine ,
Adjustable Spacing HIFU Beauty Machine ,
ABS Ultrasound Beauty Device**



Product Description

Professional Focused Ultrasound Beauty Device With Adjustable Spacing

Product Description:

The HIFU Beauty Machine boasts an impressive energy output range of 0.1J-2.5J (0.1J/step), which can be easily adjusted to suit individual client requirements. With its 300,000 shots (IPL) lamp longevity, this machine is built to last and provide consistent results. The pack size of the HIFU Beauty Machine is 60*44*46cm, making it compact and easy to store when not in use. Its spacing feature is also customizable, with a range of 1.0-10mm (0.5mm/step) that can be adjusted to suit the specific needs of each client. Investing in the HIFU Beauty Machine means investing in the latest technology for non-invasive skin tightening beauty therapy. This machine is user-friendly, efficient, and delivers results that will keep clients coming back for more.

Features:

Product Name: HIFU Beauty Machine

Energy Output: 0.1J-2.5J (0.1J/step)

Body Type: Vertical

Specification: Packing, 64 X 81 X 131 Cm


Size: Portable

Scope: 5.0mm-25mm

Keywords: Hifu Facial Machine, Home-use HIFU Beauty Device, HIFU Face Lift Device

Technical Parameters:

Skin Tightening Beauty Device	Hifu Ultrasound Facelift Machine
Lamp Longevity	300,000 Shots (IPL)
Specification	Packing, 64 X 81 X 131 Cm
Body Type	Vertical
Supplier Info	Manufacturer
Spacing	1.0-10mm (0.5mm/step)
Energy Output	0.1J-2.5J (0.1J/step)
Material	ABS
Energy Type	HIFU (High Intensity Focused Ultrasound)
Input Voltage	AC100V-240V
Working Mode	Simultaneous / Independent Control

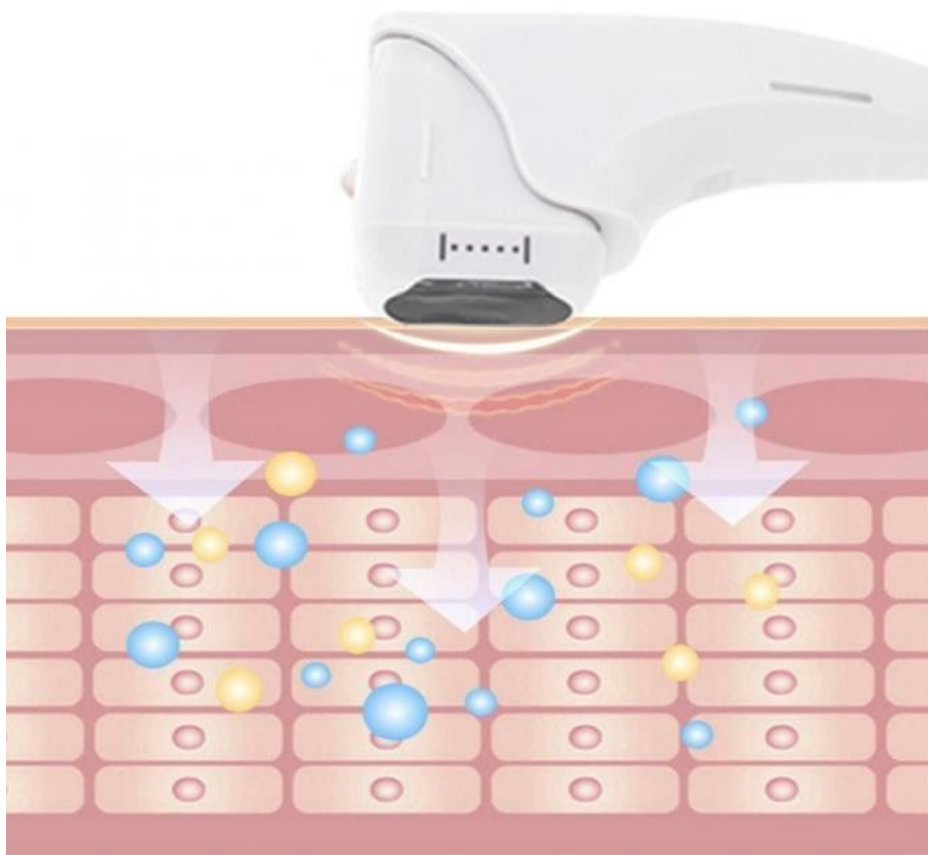

Tight, delicate and not loose

Create a small "V" face, scientific 45 elevation angle, more suitable for facial contour, say goodbye to loose skin



no fear of time Reshape youthful-looking muscles

Reasons for Choice

guard young Rejuvenate the skin

wrinkle

Anti-wrinkle/smooth fine lines

shrink

Shrink/reduce pores

Lift and tighten

Firming/enhancing elasticity

shape


Sculpting/Bumping



Rejuvenate the skin and break the "wrinkle" of the skin

Through high-frequency radio waves, the heat penetrates into the inner layer of the skin, awakening the vitality of collagen and caring for the skin.

III

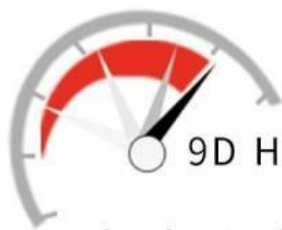

9D HIFU



1 minute / 12 rows of RBI
Multi-line dots save half the time

Up to 12 rows of points at a time



9D HIFU



Hit 12 points at a time


Ordinary HIFU



Hit 1 row at a time

Multiple scene selection

Meet various personnel needs

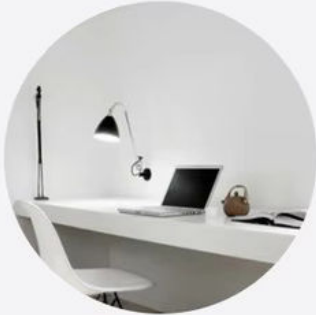

workplace

+



beauty club



personal

+



family

Customization:


If you're looking for a high-quality ultrasound beauty equipment, consider the HIFU-3 3D HIFU Machine by other. This HIFU Face Slimming Device is manufactured in Guangzhou, China and has CE certification. The minimum order quantity is 1 set with a discussible price and regular packaging. Delivery time is approximately 5-8 working days. Payment terms include L/C and T/T, and the supply ability is up to 10,000 pieces per month.

The HIFU-3 3D HIFU Machine has an energy output range of 0.1J-2.5J (0.1J/step) and is manufactured by a reputable supplier. It has two working modes, simultaneous and independent control, and is suitable for all skin types. The spacing is adjustable from 1.0-10mm (0.5mm/step).

Working Principle


- Vaginal Tightening

It uses ultrasonic waves to send out ultrasonic energy focusing in the lamina propriety and muscle fiber layer in a predetermined depth. In 0.1 second, the temperature of the region can reach to above 65°C, so the collagen is reorganized and the normal tissue outside the focal region is undamaged. The desired depth layer can obtain the ideal effect of collagen contraction, reorganization and regeneration.

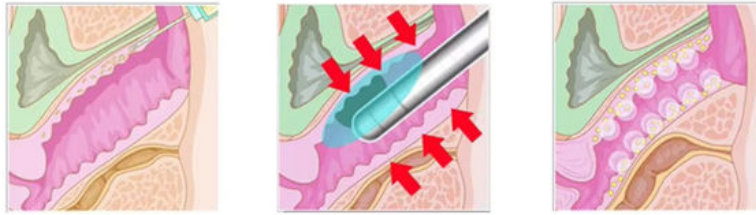

- Facail lifting

HIFU technology will release a highly concentrated focus of ultrasonic energy and deliver into the deep skin and cause thermal coagulation to SMAS layer. High heat is condensed in the appropriate location, to stimulate the skin deep dermal tissue, produce more collagen, so it make the skin more compact.

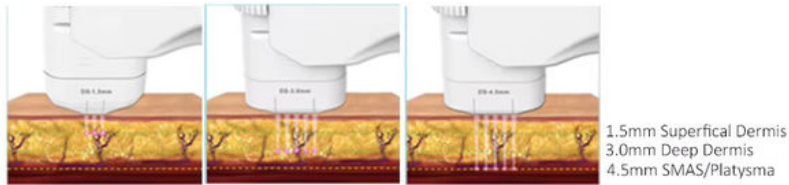

OUR EQUIPMENT



Product parameters

Model	HU700
Handle	3 Handles
Screen	15inch color touch screen
Energy Source	High Intensity Focused Ultrasound
Shot Number (3D HIFU Handle)	20 000 shots/each cartridge (1.5mm/3.0mm/4.5mm/6.0mm/8.0mm/10mm/13mm/16mm)
Shot Number (4D HIFU Handle)	62 000 shots/each cartridge (1.5mm/3.0mm/4.5mm)
Shot Number (Vaginal HIFU Handle)	10 000 shots/each cartridge (3.0mm/4.5mm)
Energy for face/body	0.2 - 2 J(0.1J/step)
Cartridge	Face & Body cartridge: 10 MHz -- 1.5mm 7 MHz -- 3.0mm / 4.5mm 4 MHz -- 6.0mm / 8.0mm / 10mm / 13mm / 16mm
Energy for vaginal	0.2-2.0J
Size	machine size (35.5cm*30.5cm*23cm), package size (70cm*58cm*60cm)
Net Weight	9kg

Support and Services:

The HIFU Beauty Machine is a cutting-edge device that uses high-intensity focused ultrasound technology to provide non-invasive, safe, and effective skin rejuvenation and lifting treatments. Our product technical support and services include:

- In-depth product training and user manuals to ensure safe and effective use of the device
- One-year warranty to protect your investment
- Ongoing technical support to answer any questions or concerns you may have
- Regular software updates to ensure optimal performance and results
- Marketing and promotional materials to help you promote and grow your business

Our team of experts is dedicated to providing the highest level of support and service to ensure your success with the HIFU Beauty Machine.

Packing and Shipping:

Product Packaging:

- 1 HIFU Beauty Machine
- 1 Power Adapter
- 1 User Manual

Shipping:

- We offer free shipping worldwide.
- Orders are processed within 24 hours and shipped within 1-3 business days.
- Delivery times vary depending on location, but typically range from 7-14 business days.

FAQ:

Q1: What is the brand name of the HIFU Beauty Machine?

A1: The brand name of this machine is other.

Q2: What is the model number of the HIFU Beauty Machine?

A2: The model number of this machine is hifu-3.

Q3: Where is the HIFU Beauty Machine manufactured?

A3: The HIFU Beauty Machine is manufactured in Guangzhou, China.

Q4: Is the HIFU Beauty Machine certified?

A4: Yes, this machine is CE certified.

Q5: What is the minimum order quantity for the HIFU Beauty Machine?

A5: The minimum order quantity for this machine is 1 set.

Q6: What is the price for the HIFU Beauty Machine?

A6: The price for this machine is discussible and can vary based on quantity and other factors.

Q7: What are the packaging details for the HIFU Beauty Machine?

A7: The HIFU Beauty Machine comes with regular packaging.

Q8: How long does delivery of the HIFU Beauty Machine take?

A8: Delivery time for this machine is 5-8 working days.

Q9: What payment terms are accepted for the HIFU Beauty Machine?

A9: The accepted payment terms are L/C and T/T.

Q10: What is the supply ability of the HIFU Beauty Machine?

A10: The supply ability of this machine is 10000 pieces per month.



Guangzhou Lizhirun Technology Development Co., Ltd



18688905627



cenglinlin25@gmail.com



hifubeaity-machine.com

Room 502, 5th Floor, Building F, Huixin Industrial Zone, No. 89, Hetai Road, Hebain Village, Baiyun District,
Guangzhou